

Operating Buk Bedienung Bediening Betjening Emploi




Cleaning and maintenance


Stainless steel surfaces should be cleaned with a microfibre cloth, glass cleaner or other non abrasive and gentle agents.

Attention:

Never use solvents or aggressive household cleaners like bleach, acids (e.g. descalers) or ammonia-based bathroom cleaner to clean the product. Steel wool and hard brushes can permanently damage the surfaceof the steel.

COVER FOR ARM SUPPORT - N1961-1 (Additional use)

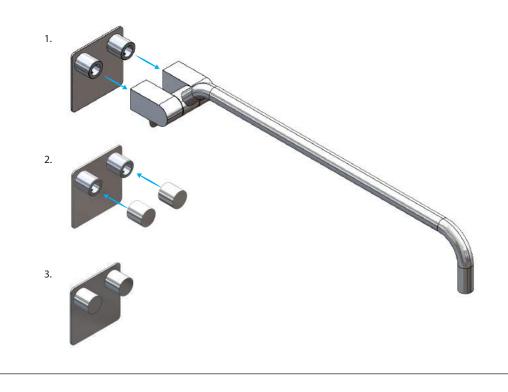

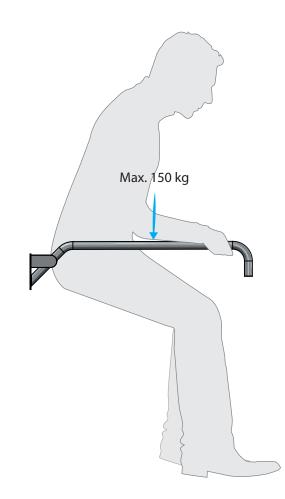

Tested by everyday life

To ensure that the arm rest is capable of withstanding heavy load and conditions every day, the product has been tested in accordance to the requirements of the European standard for technical aids for disabled, BS/EN 12182:2012. This includes a static proof load test with 150 kg, following by a static overload test with 225 kg.

We recommend all FROST CARE products are installed on a solid wall. Reinforcement should always be used if the wall is not stable. All tests are carried out on products mounted on a solid wall, thus we cannot guarantee the same performance and durability if this is not observed.

This product has a maximum weight capacity of 150 kg.

For safety reasons, we recommend to check all bolt connections to ensure that no parts are loose and that all attachment to the wall are intact.

